



## Nori® Guinea Pig Endostatin ELISA Kit- DataSheet

Endostatin is a naturally-occurring, 20-kDa C-terminal fragment derived from type XVIII collagen and is a broad-spectrum angiogenesis inhibitor and may interfere with the pro-angiogenic action of growth factors such as bFGF/FGF-2 and VEGF.<sup>[1]</sup> It was first found secreted in the media of non-metastasizing mouse cells from a hemangioendothelioma cell line in 1997 and was subsequently found in humans.<sup>[2][3]</sup> It is produced by proteolytic cleavage of collagen XVIII, a member of the multiplexin family that is characterized by interruptions in the triple helix creating multiple domains, by proteases such as cathepsins.<sup>[4]</sup> Endostatin, as a fragment of collagen 18, demonstrates a role of the ECM in suppression of neoangiogenesis.<sup>[2]</sup> Overall, endostatin down regulates many signaling cascades like ephrin, TNF- $\alpha$ , and NF $\kappa$ B signaling as well as coagulation and adhesion cascades.<sup>[5]</sup> In-vitro studies have shown endostatin blocks the proliferation and organization of endothelial cells into new blood vessels.<sup>[6]</sup> In animal studies endostatin inhibited angiogenesis and growth of both primary tumors and secondary metastasis.<sup>[2]</sup> Endostatin is currently being studied as part of cancer research. Prior results indicated that endostatin can be beneficial in combinations with other medicines, but endostatin alone gave no significant improvements in tumor/disease progression. Endostatin may also be useful as a therapeutic for inflammatory diseases like rheumatoid arthritis as well as Crohn's disease, diabetic retinopathy, psoriasis, and endometriosis by reducing the infiltration of inflammatory cells through invading angiogenesis.<sup>[7]</sup> Down's syndrome patients seem to be protected from diabetic retinopathy due to an additional copy of chromosome 21, and elevated expression of endostatin.<sup>[8]</sup>

### References

1. Folkman, J. (2006). *Exp. Cell. Res.* **312** (5, part 2): 594–607.
2. O'Reilly, M.S.; et al. (1997). *Cell* **88** (2): 277–85.
3. Standker, L; et al. (1997). *FEBS Lett* **420** (2–3): 129–33.
4. Felbor, U; et al. (2000). *EMBO J* **19** (6): 1187–94.
5. Abdollahi, A.; et al. (2004). *Mol. Cell* **13** (5): 649–63.
6. Folkman, J.; Kalluri, R. (2004). "Cancer without disease". *Nature* **427** (6977): 787.
7. Yin, G., et al. (2002). *Mol Ther* **5** (5 Pt 1): 547–54.
8. Ryeom S, (2009). *J Craniofac Surg.* 20 Suppl 1: 595–6.

### PRINCIPLE OF THE ASSAY

This ELISA kit is for quantification of endostatin in Guinea Pig samples. This is a quick ELISA assay that reduces time to 50% compared to the conventional method, and the entire assay only takes 3 hours. This assay employs the quantitative sandwich enzyme immunoassay technique and uses biotin-streptavidin chemistry to improve the performance of the assays. An antibody specific for Guinea Pig endostatin has been pre-coated onto a microplate. Standards and samples are pipetted into the wells and any endostatin present is bound by the immobilized antibody. After washing away any unbound substances, a detection antibody specific for Guinea Pig endostatin is added to the wells. Following wash to remove any unbound antibody reagent, a detection reagent is added. After intensive wash a substrate solution is added to the wells and color develops in proportion to the amount of endostatin bound in the initial step. The color development is stopped and the intensity of the color is measured.

This package insert must be read in its entirety before using this product.

### Storage

Store at 4 °C. The kit should be used in 12 months.



## Nori® Guinea Pig Endostatin ELISA Kit- DataSheet

### MATERIALS PROVIDED

Description	Quantity	Description	Quantity	Description	Quantity
Antibody Precoated Plate	1	20 x PBS	1	Substrate Solution	1
Detection Antibody	1	20 x Assay Buffer	1	Stop Solution	1
HRP Conjugate	1	Reagent Diluent	1	DataSheet/Manual	1
Standard	3	MSDS/CoA	1	96-well plate sheet	1

Bring all reagents to room temperature before use.

### Reagent Preparations

**Guinea Pig Endostatin Detection Antibody** (1 vial) – The lyophilized Detection Antibody should be stored at -20°C in a manual defrost freezer for up to 12 months, if not used immediately. Centrifuge for 1 min at 6000 x g to bring down the material prior to open the vial. The vial contains sufficient Detection Antibody for a 96-well plate. Add 200 µL of sterile 1 x PBS to a vial and vortex 30 and allow it to sit for 5 min prior to use. Make 1:50 dilution in Reagent Diluent. Take 200 µL of detection antibody to 10.5 mL of Reagent Diluent to make **Working dilution of detection antibody** and vortex 30 sec if the entire 96-well plate is used. If the partial antibody is used store the rest at -20°C until use.

**Guinea Pig Endostatin Standard** (3 vials) – The lyophilized Guinea Pig Endostatin Standard has a total of 3 vials. Each vial contains the standard sufficient for generating a standard curve. The unreconstituted standard can be stored at -20°C for up to 12 months if not used immediately. Centrifuge for 1 min at 6000 x g to bring down the material prior to open the vial. Add 500 µL of 1 x Assay Buffer to a vial to make the high standard concentration of 8,000 pg/ml and **vortex 1 min** and allow it to sit for 5 min. A seven-point standard curve is generated using 2-fold serial dilutions in 1 x Assay Buffer, **vortex 30 sec** for each of dilution step.

**HRP Conjugate** (53 µL) – Centrifuge for 1 min at 6000 x g to bring down the material prior to open the vial. The vial contains 53 µL HRP Conjugate sufficient for one 96-well plate. If the volume is less than 53 µL, add sterile 1 x PBS to reach 53 µL and vortex 10 sec. Make 1:200 dilutions in Reagent Diluent. If the entire 96-well plate is used, add 53 µL of HRP Conjugate to 10.5 mL of Reagent Diluent to make **working dilution of HRP Conjugate** and vortex 30 sec prior to the assay. The rest of undiluted HRP Conjugate can be stored at 4°C for up to 12 months. DO NOT FREEZE.

**20 x PBS**, pH 7.3, 25 mL- Dilute to 1 x PBS with deionized distilled water and mix well prior to use.

**20 x Assay Buffer**, 20 mL- Dilute to 1 x Assay Buffer with 1 x PBS prior to use.

**Reagent Diluent**, 21 mL.

**Substrate Solution**, 10.5 mL.

**Stop Solution**, 5.5 mL.



## Nori® Guinea Pig Endostatin ELISA Kit- DataSheet

### Assay Procedure

1. Lift the plate cover from the top left and cover the wells that are not used. **Vortex the standards and samples for 10 sec** before applying to the plate. Add 100 µL of **diluted sample** (see below) or **standard** per well and use duplicate wells for each standard or sample. Cover the 96-well plate and incubate 1 hour at room temperature. Attention: **MUST vortex standards and samples for 10 sec before pipetting to the wells!**
2. Aspirate each well and wash with **1 x Assay Buffer**, repeating the process one time for a total of two washes. Wash by filling each well with 1 x Assay Buffer (300 µL) using a multi-channel pipette, manifold dispenser or auto-washer. Complete removal of liquid at each step is essential for good performance. After the last wash, remove any remaining Assay Buffer by aspirating or by inverting the plate and blotting it against clean paper towels.
3. Add 100 µL of the **working dilution of Detection Antibody** to each well. Cover the plate and incubate 1 hour at room temperature.
4. Repeat the aspiration/wash as in step 2.
5. Add 100 µL of the **working dilution of HRP Conjugate** to each well. Cover the plate and incubate for 20 minutes at room temperature. Avoid placing the plate in direct light.
6. Repeat the aspiration/wash as in step 2.
7. Add 100 µL of **Substrate Solution** to each well. Incubate for up to 20 minutes (**depending on signal**) at room. Over incubation of the substrate will result in over flow of high standard and thus should be avoided. Avoid placing the plate in direct light.
8. Add 50 µL of **Stop Solution** to each well. Gently tap the plate to ensure thorough mixing.
9. Determine the optical density of each well immediately, using a microplate reader set to 450 nm. If wavelength correction is available, set to 540 nm or 570 nm. If wavelength correction is not available, subtract readings at 540 nm or 570 nm from the readings at 450 nm. This subtraction will correct for optical imperfections in the plate. Readings made directly at 450 nm without correction may be higher and less accurate.

**Sample dilution:** Samples should be diluted with four volumes of 1 x Assay Buffer and **vortex for 1 min** prior to assay. If the OD value still exceeds the upper limit of the standard curve, further dilution is recommended till it falls in the detection range and the dilution factor must be used for calculation of the concentration.

### **Precaution and Technical Notes**

1. It is critical to follow the procedure step by step otherwise appropriate color development may not occur as expected.
2. A standard curve should be generated for each set of samples assayed. Thorough mixing of standards at each of dilution steps is critical to acquire a normal standard curve and **vortex again (10 sec) before pipetting to the 96-well plate.**
3. HRP Conjugate contains enzyme, DO NOT mass up with Detection Antibody.
4. The Stop Solution is an acid solution, handle with caution.
5. This kit should not be used beyond the expiration date on the label.
6. A thorough and consistent wash technique is essential for proper assay performance.
7. Use a fresh reagent reservoir and pipette tips for each step.
8. It is recommended that all standards and samples be assayed in duplicate.
9. Avoid microbial contamination of reagents and buffers. This may interfere with the sensitivity of the assay.

## Nori® Guinea Pig Endostatin ELISA Kit- DataSheet


### Calculation of Results

Average the duplicate readings for each standard, control, and sample and subtract the average zero (blank) standard optical density.

Create a standard curve by reducing the data using computer software capable of generating a four-parameter logistic (4-PL) curve-fit. As an alternative, construct a standard curve by plotting the mean absorbance for each standard on the y-axis against the concentration on the x-axis and draw a best fit curve through the points on the graph. The data may be linearized by plotting the log of the endostatin concentrations versus the log of the O.D. and the best fit line can be determined by regression analysis. This procedure will produce an adequate but less precise fit of the data. If samples have been diluted, the concentration read from the standard curve must be multiplied by the dilution factor.

### The Standard Curve

The graph below represents typical data generated when using this Guinea Pig endostatin ELISA Kit. The standard curve was calculated using a computer generated 4-PL curve-fit. For this case, a Bio-Rad iMark™ Microplate Reader and a Microplate Manager 6 Software were used to generate this curve. The correlation coefficient ( $r^2$ ) is 0.999-1.000.

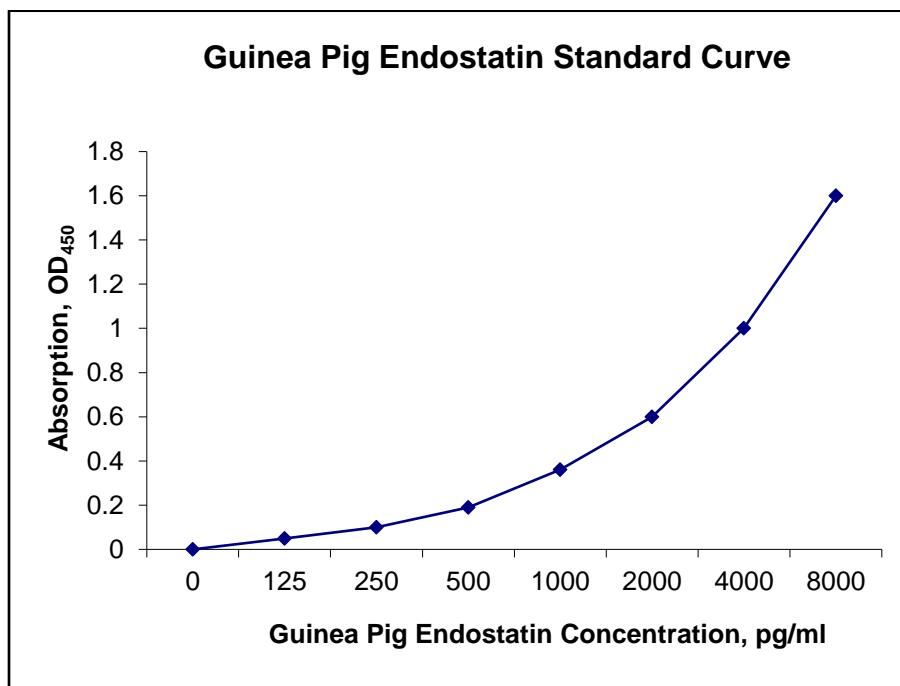



## Nori® Guinea Pig Endostatin ELISA Kit - DataSheet

### Specificity

The following recombinant Guinea Pig proteins prepared at 1 ng/ml were tested and exhibited no cross-reactivity or interference.

Adiponectin, ApoAI, BMP1, BMP2, BMP3, BMP4, BMP5, BMP6, BMP7, CCL2, CCL4, CCL5, CRP, FGF acidic, HGF, HSP27, IGF-1, IL-1 $\alpha$ , IL-1 $\beta$ , IL-1ra, IL-4, IL-5, IL-6, IL-8, IL-10, IL-12, IL-13, IL-15, IL-17C, ENDOSTATIN1, ENDOSTATIN3, IFN- $\alpha$ , IFN- $\beta$ , IFN $\gamma$ , MMP-2, MMP-3, MMP-9, PDGF, PLA2G7, prolactin, TLR1, TLR2, TLR3, TLR4, TLR9, TGF- $\beta$ 1, TGF- $\beta$ 2, TGF- $\beta$ 3, TNF- $\alpha$ , TNF RI, TNF RII, VEGF, VEGF-R1.

### Calibration

This kit is calibrated against a highly-purified yeast-expressed recombinant Guinea Pig endostatin.

### Detection Range

125-8000 pg/ml

### Assay Sensitivity

25 pg/ml

### Assay Precision

Intra-Assay %CV: 5; Inter-Assay %CV: 10

**Sample Types:** Plasma, serum, cell/tissue lysates, cell culture supernatant, synovial fluid (SF), bronchoalveolar lavage (BAF), cerebrospinal fluid (CSF), urine, and other biological fluid.

## For Research Use Only

### Related products

20 x sample diluent, GR103061

10 x ELISA Wash Buffer, GR103014

10 x Reagent Diluent, GR103028

20 x PBS, GR103004

ELISA Substrate, GR103021

ELISA Stop Solution, GR103055

ELISA HRP Conjugate, GR103044

Guinea Pig endostatin standard

Guinea Pig endostatin detection antibody



## Nori® Guinea Pig Endostatin ELISA Kit DataSheet

### Troubleshooting Guide

Problem	Possible causes	Solution
Poor standard curve	<ul style="list-style-type: none"><li>• Inaccurate pipetting</li><li>• Improper standard curve</li><li>• OD<sub>450</sub> too high for the high standard point</li></ul>	<ul style="list-style-type: none"><li>• Check pipettes</li><li>• Check and use the correct dilution buffer</li><li>• Vortex 30 sec for each of standard dilution steps and <b>vortex again (10 sec) before pipetting to the 96-well plate.</b></li><li>• Reduce substrate incubation time</li></ul>
Low signal	<ul style="list-style-type: none"><li>• Improper preparation of standard, samples, detection antibody, and/or HRP Conjugate</li><li>• Too brief incubation times</li><li>• Inadequate reagent volume or improper dilution</li></ul>	<ul style="list-style-type: none"><li>• Briefly spin down vials before opening. Reconstitute the powder thoroughly.</li><li>• Ensure sufficient incubation time including substrate incubation.</li><li>• Check pipettes and ensure correct preparation.</li></ul>
Large CV	<ul style="list-style-type: none"><li>• Inaccurate pipetting and mixing</li><li>• Improper standard/sample dilutions.</li><li>• Air bubbles in wells.</li></ul>	<ul style="list-style-type: none"><li>• Check pipettes and ensure thorough mixing.</li><li>• Use the correct dilution buffers</li><li>• Remove bubbles in wells.</li></ul>
High background	<ul style="list-style-type: none"><li>• Plate is insufficiently washed.</li><li>• Contaminated wash buffer</li></ul>	<ul style="list-style-type: none"><li>• Review the datasheet for proper wash. If using a plate washer, ensure that all ports are unobstructed.</li><li>• Make fresh wash buffer</li></ul>
No signal detected	<ul style="list-style-type: none"><li>• The procedure was misconducted.</li></ul>	<ul style="list-style-type: none"><li>• Ensure the step-by-step protocol was correctly followed and no misstep was conducted.</li></ul>
Low sensitivity	<ul style="list-style-type: none"><li>• Improper storage of the ELISA kit</li><li>• Stop solution</li></ul>	<ul style="list-style-type: none"><li>• Store standards and detection antibody at -20°C after reconstitution, others at 4°C. Keep substrate protected from light.</li><li>• Adding stop solution to each well before reading plate</li></ul>