

Nori® Equine NT-proBNP ELISA Kit DataSheet

Brain natriuretic peptide is a 32-amino acid polypeptide secreted by the ventricles of the heart in response to excessive stretching of cardiomyocytes. BNP is synthesized as a 134-amino acid preprohormone (preproBNP), encoded by the gene NPPB.[1] Removal of the 25-residue N-terminal signal peptide generates the prohormone, proBNP, which is stored intracellularly as an O-linked glycoprotein; proBNP is subsequently cleaved between arginine-102 and serine-103 by a specific convertase into NT-proBNP and the biologically active 32-amino acid polypeptide BNP-32, which are secreted into the blood in equimolar amounts.[2] The release of BNP is modulated by calcium ions. BNP is secreted attached to a 76-amino acid N-terminal fragment in the prohormone called NT-proBNP (BNPT), which is biologically inactive. Once released, BNP binds to and activates the atrial natriuretic factor receptors NPRA, and to a lesser extent NPRB, in a fashion similar to atrial natriuretic peptide (ANP) but with 10-fold lower affinity. The biological half-life of BNP, however, is twice as long as that of ANP, and that of NT-proBNP is even longer, making these peptides better targets than ANP for diagnostic blood testing. The physiologic actions of BNP are similar to those of ANP and include decrease in systemic vascular resistance and central venous pressure as well as an increase in natriuresis. The net effect of these peptides is a decrease in blood pressure due to the decrease in systemic vascular resistance and, thus, afterload. Additionally, the actions of both BNP and ANP result in a decrease in cardiac output due to an overall decrease in central venous pressure and preload as a result of the reduction in blood volume that follows natriuresis and diuresis. A normal level BNP or NT-proBNP can rule out acute heart failure in the emergency setting. However, an elevated BNP or NT-proBNP should never be used to "rule in" acute or chronic heart failure in the emergency setting due to lack of specificity.[3] BNP and NT-proBNP are also typically increased in patients with left ventricular dysfunction, with or without symptoms (BNP accurately reflects current ventricular status, as its half-life is 20 minutes, as opposed to 1–2 hours for NT-proBNP).[4] Low BNP was found to be a predictor of survival to age 90 in men,[5] and BNP may be a reliable predictor of cardiovascular mortality in diabetics.[18]

References

1. Sudoh T, et al. (1989) *Biochem. Biophys. Res. Commun.* 159 (3), 1427-1434.
2. Schellenberger U, et al. (2006). *Arch. Biochem. Biophys.* 451 (2): 160–6.
3. Maisel A, et al. (2002) *N Engl J Med.* 347 (3): 161–7.
4. Atisha D, et al. (2004) *Am. Heart J.* 148 (3): 518–23.
5. Nilsson, G; e tal. (2014) *Healthy Aging Research.* 3 (5): 1–10.
6. Bhalla MA, et al. (2004). *J. Am. Coll. Cardiol.* 44 (5): 1047–52.

PRINCIPLE OF THE ASSAY

This is a quick ELISA assay that reduces time to 50% compared to the conventional method, and the entire assay only takes 3 hours. This assay employs the quantitative sandwich enzyme immunoassay technique and uses biotin-streptavidin chemistry to improve the performance of the assays. An antibody specific for NT-proBNP has been pre-coated onto a microplate. Standards and samples are pipetted into the wells and any NT-proBNP present is bound by the immobilized antibody. After washing away any unbound substances, a detection antibody specific for NT-proBNP is added to the wells. Following wash to remove any unbound antibody reagent, a detection reagent is added. After intensive wash a substrate solution is added to the wells and color develops in proportion to the amount of NT-proBNP bound in the initial step. The color development is stopped, and the intensity of the color is measured.

This package insert must be read in its entirety before using this product.

Storage Store at 4 °C. The kit can be used in 6 months.



Nori® Equine NT-proBNP ELISA Kit DataSheet

MATERIALS PROVIDED

Description	Quantity	Description	Quantity	Description	Quantity
Antibody Precoated Plate	1	20 x PBS	1	Substrate Solution	1
Detection Antibody	1	20 x Assay buffer	1	Stop Solution	1
HRP Conjugate	1	96-well plate sheet	1	DataSheet	1
Standard	3				

Bring all reagents to room temperature before use.

Reagent Preparations

Equine NT-proBNP Detection Antibody (1 vial) – The lyophilized Detection Antibody should be stored at 4°C in a manual defrost freezer for up to 6 months, if not used immediately. Centrifuge for 1 min at 6000 x g to bring down the material prior to open the vial. The vial contains sufficient Detection Antibody for a 96-well plate. Add 200 µL of sterile 1 x PBS to the antibody vial and vortex briefly and sit for 5 min. Take 200 µL of detection antibody to 10.5 mL of 1 x PBS to make **working dilution of Detection Antibody** if the entire 96-well plate is used. If the partial antibody is used store the rest at - 20°C until use.

Equine NT-proBNP Standard (3 vials) – Each lyophilized Equine NT-proBNP Standard vial contains the standard sufficient for generating a standard curve. The non-reconstituted standard can be stored at 4°C for up to 6 months if not used immediately. Centrifuge for 1 min at 6000 x g to bring down the material prior to open the tube. Add 500 µL of 1 x Assay Buffer to one Standard vial to make the high standard concentration of 50 ng/mL. Vortex 20 sec and allow it to sit for 5 min prior to use. A seven-point standard curve is generated using 2-fold serial dilutions in the Assay Buffer, vortex 20 sec for each of dilution step.

HRP Conjugate (55 µl) – Centrifuge for 1 min at 3000 x g to bring down the material prior to open the vial. The vial contains 55 µL HRP Conjugate sufficient for one 96-well plate (Please notify us if it is below 40 µL). Add 1 x PBS to reach 55 µL if needed. Make 1:200 dilutions in 1 x PBS. If the entire 96-well plate is used, add all HRP Conjugate to 10.5 mL of 1 x PBS to make **working dilution of HRP Conjugate** and vortex 30 sec prior to the assay. The rest of undiluted HRP Conjugate can be stored at 4 °C for up to 6 months. It is recommended to fully recover the HRP with 1 ml of 1 x PBS at final use.

20 x PBS, pH 7.3, 25 mL- Dilute to 1 x PBS with deionized distilled water and mix well prior to use.

20 x Assay buffer, 20 mL- Dilute to 1 x Assay buffer with 1 x PBS prior to use.

Substrate Solution, 10.5 mL.

Stop Solution, 5.5 mL.

Nori® Equine NT-proBNP ELISA Kit DataSheet

1. All procedures are conducted at room temperature (20-25 °C) and ensure **equal pipetting/dispensing** at each step and remove air bubbles in the wells for all steps.
2. Lift the plate cover and cover the unused wells or reseal the unused strips in the aluminum bag with desiccant at 4 °C. **Vortex the standards and samples for 10 sec** before applying to the plate. Add 100 µL of **sample** or **standard** per well and use duplicate wells for each standard or sample. Cover the 96-well plate and incubate on shaker at 450 rpm for **1 h** (1.5 h if no shaking). Attention: **MUST vortex standards and samples for 10 sec before pipetting to the wells!**
3. Aspirate each well and wash with 300 µL of **1 x Assay Buffer** for two times. Wash by filling each well with 1 x Assay Buffer using a multi-channel pipette, manifold dispenser, squirt bottle or auto-washer. Complete removal of liquid at each step is essential for good performance. After the last wash, remove any remaining Assay Buffer by aspirating or by inverting the plate and blotting it against clean paper towels.
4. Add 100 µL of the **working dilution of Detection Antibody** to each well. Cover the plate and incubate on shaker at 450rpm for 1 h (1.5 h if no shaking).
5. Repeat the aspiration/wash as in step 3.
6. Add 100 µL of the **working dilution of HRP Conjugate** to each well. Cover the plate and incubate for 20 minutes. Avoid placing the plate in direct light.
7. Repeat the aspiration/wash as in step 3 but **wash 4 times** instead.
8. Add 100 µL of **Substrate Solution** to each well and observe the color development every 1-2 mins. Incubate for up to **30 minutes depending on signal** (**Stop** the reaction when it turns to dark blue in the highest standard well). Over-incubation of the substrate will result in overflow of high standard and should be avoided. Avoid placing the plate in direct light.
9. When it gets to dark blue in the highest concentration of standard wells, add 50 µL of **Stop Solution** to each well to stop the reaction. Gently tap the plate to ensure thorough mixing. **Ensure all wells turn yellow by pipette tip prior to measurement.**
10. Determine the optical density of each well immediately, using a validated microplate reader set to 450 nm. If wavelength correction is available, set to 540 nm or 570 nm. If wavelength correction is not available, subtract readings at 540 nm or 570 nm from the readings at 450 nm. This subtraction will correct for optical imperfections in the plate. Readings made directly at 450 nm without correction may be higher and less accurate.

Sample dilution: If high density is expected, samples should be diluted with equal volume of 1 x Assay Buffer and **vortex for 1 min** prior to assay. If the OD value still exceeds the upper limit of the standard curve, further dilution is recommended till it falls in the detection range and the dilution factor must be used for calculation of the concentration.

Precaution and Technical Notes

1. It is critical to follow the procedure step by step otherwise appropriate color development may not occur as expected and make sure no air bubbles in wells before adding reagents.
2. A standard curve should be generated for each set of samples assayed. Thorough mixing of standards at each of dilution steps is critical to acquire a normal standard curve and **vortex again (10 sec) before pipetting to the 96-well plate.**
3. HRP Conjugate contains enzyme, DO NOT mass up with Detection Antibody.
4. The Stop Solution is an acid solution, handle with caution.
5. This kit should not be used beyond the expiration date on the label.
6. A thorough and consistent wash technique is essential for proper assay performance.
7. Use a fresh reagent reservoir and pipette tips for each step.
8. It is recommended that all standards and samples be assayed in duplicate.
9. Avoid microbial contamination of reagents and buffers. This may interfere with the performance of the assay.

Nori® Equine NT-proBNP ELISA Kit DataSheet


Calculation of Results

Average the duplicate readings for each standard, control, and sample and subtract the average zero (blank) standard optical density.

Create a standard curve by reducing the data using computer software capable of generating a four-parameter logistic (4-PL) curve-fit. As an alternative, construct a standard curve by plotting the mean absorbance for each standard on the y-axis against the concentration on the x-axis and draw a best fit curve through the points on the graph. The data may be linearized by plotting the log of the NT-proBNP concentrations versus the log of the O.D. and the best fit line can be determined by regression analysis. This procedure will produce an adequate but less precise fit of the data. If samples have been diluted, the concentration read from the standard curve must be multiplied by the dilution factor.

The Standard Curve

The graph below represents typical data generated when using the ELISA Kit. The correlation coefficient (r^2) is 0.999-1.000. The standard curve should be calculated using a computer generated 4-PL curve-fit to determine concentrations of unknown specimens.

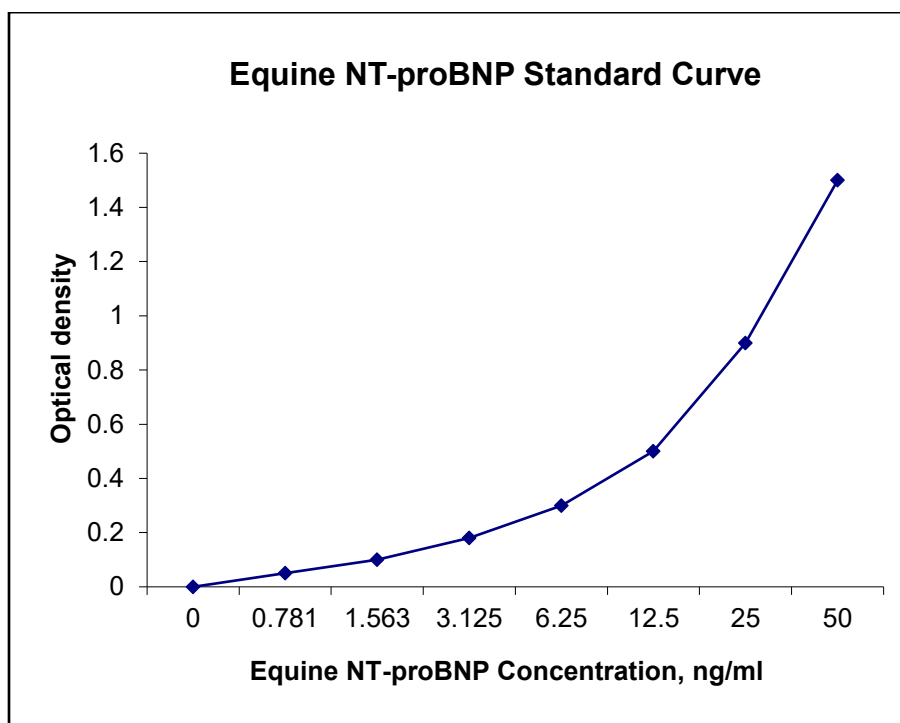



Nori® Equine NT-proBNP ELISA Kit DataSheet

Specificity

The following recombinant Equine proteins prepared at 10 ng/ml were tested and exhibited no cross-reactivity or interference.

ApoA1, BMP4, CCL4/MIP-1 β , CRP, HSP27, IL-1 β , IFN γ , PDGF, PLA2G7, prolactin, TGF β 1, TLR3, TNF- α , TNF RI, TNF RII, VEGF.

Detection Range

0.78-50 ng/ml

Assay Sensitivity

150 pg/ml

Assay Precision

Intra-Assay %CV: 5; Inter-Assay %CV: 8

Related products

1. GR238016 50 ml Reagent Reservoir, 100/case, 5 packs/case (pack of 20)
2. GR238004 Tissue Culture 96-well Microplate, individually packed, Case of 50
3. GR238019 1.5 ml Microcentrifuge tube with screw cap and free-standing, pack of 500
4. GR238007 125 ml leak-resistant HDPE bottle, colorless, pack of 24
5. GR238002 Microplate 12x8-Well Strip High Binding, Case of 50
6. GR238003 Microplate 12x8-Well Strip Medium Binding, Case of 50
7. GR238032 42592 Costar Stripwell Microplate 1 x 8 Flat Bottom, High Binding, Case of 100
8. GR238001 468667 Thermo Microplate 12x8-Well Strip Nunc Maxisorp F8, Case of 60
9. GR238031 96-well microplate sealer plastic, pack of 100

DECLARATION

THIS REAGENT IS FOR IN VITRO LABORATORY TESTING AND RESEARCH USE ONLY. DO NOT USE IT FOR CLINICAL DIAGNOSTICS. DO NOT USE OR INJECT IT IN HUMANS AND ANIMALS.

**FOR LABORATORY RESEARCH USE ONLY
NOT FOR USE IN HUMANS AND ANIMALS**

Nori® Equine NT-proBNP ELISA Kit DataSheet

Troubleshooting Guide

Problem	Possible causes	Solution
Poor standard curve	- Inaccurate pipetting - Insufficient vortexing - OD₄₅₀ too high for the high standard point - Air bubbles in wells. - Standard defect or not fully recovered	- Check pipette calibration and ensure equal dispensing. - Vortex 30 sec for each of standard dilution steps and vortex again (10 sec) before pipetting to the 96-well plate. - Reduce substrate incubation time - Remove air bubbles in wells by pipette tip. - Change a standard vial or spin down the vial before reconstitution
Low signal	- Improper preparation of reagents and storage - Too brief incubation times - Inadequate reagent volume or improper dilution - Standard defect and sample overdiluted	- Spin down vials before opening. Reconstitute the content thoroughly. Proper storage of plate and strip and detection antibody after first usage. - Microplate shaker may improve signals. - Insufficient HRP Conjugate. Ensure sufficient incubation time and increase sample incubation to 2 h. - Change a Standard vial. Undilute sample or less dilution
Overflow in the standards	- Substrate incubation too long - Air bubbles in wells	- Observe the color development every 1-2 mins and reduce substrate incubation time. - Stop the reaction by adding 50 µl of Stop Solution when it turns to dark blue in the highest concentration of standard wells. - Remove air bubbles in wells
Large CV	- Inaccurate pipetting and mixing - Improper sample dilutions. - Air bubbles in wells. - Microplate reader out of calibration - It did not turn yellow after adding Stop Solution	- Check pipettes and ensure the pipette is calibrated properly. - Ensure accurate pipetting and thorough mixing. - Use reverse, instead of forward pipetting. - Use the correct dilution buffers - Remove air bubbles in wells by pipette tip. - Calibrate the microplate reader properly - If it did not turn yellow after adding Stop Solution, mix with pipette tip till it turns yellow prior to measurement.
High background	- Reagent reservoir issue - Plate is insufficiently washed and air bubbles in wells. - Contaminated Assay Buffer - Pipet tip contaminated	- Use a new reagent reservoir for Substrate Solution. - Increase wash to 4 times before adding substrate and ensure plate washer functions normally. Remove air bubbles in wells by pipette tip. Use squirt bottle for washing. - Make fresh Assay Buffer and wash thoroughly. - Use new pipette tips for blank wells.
No signal detected	- The procedure misconducted. - Failures of spin down the contents in Detection Antibody and Standards. - Failure of Substrate or HRP - Samples overdiluted	- Ensure the step-by-step protocol. Spin vials of Detection antibody and Standard to completely recover the content. - Ensure HRP volume. Mix 100 µl of Substrate with 0.5 µl HRP and dark blue color should develop in 5 min. - Try a new standard vial and use positive control. - Try not dilute samples
Low sensitivity	- Improper dilutions of standards - Improper storage of the ELISA kit	- Ensure accurate and thorough dilutions of standards at each step. - Store detection antibody at -20°C after reconstitution and others at 4°C. Keep substrate solution protected from light.